&me.

INSTRUCTION MANUAL MASSAGE GUN

11088811

INTRODUKTION

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Vi tager forbehold for trykfejl.

INTRODUCTION

Thank you for choosing a product from &me. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- The product engine/electronics, cord or plug cannot be submerged in water, and no water may get in contact with the engine/electronics.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept out of reach of children.

- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity, physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.

PRODUCT DISPOSAL

Used product should not be disposed as household waste. According to legislation, the product should be disposed of properly so that the materials can be recycled in an environmentally

friendly manner. Further information is available from the responsible local authority or local disposal company. The product is equipped with the following logo:

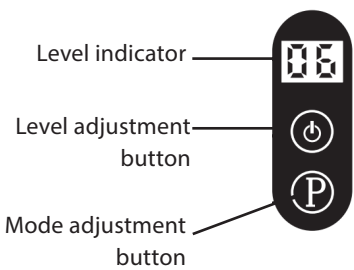
Danger: To reduce the risk of electric shock, fire, and personal injury, this product must be used according to the following instructions:

- For adults only. Use only on dry, clean skin surfaces. Do not use the device through clothing.
- Gently press and move the device on the skin. Limit use on each area to about 60 seconds.
- Use this device only on soft tissues of the body without causing pain or discomfort. Do not use it on the head or on any hard or bony parts of the body.
- Always use the vibrating head that provides the best results. Bruising may occur with moderate or normal pressure. Always check the area being treated.
- Stop using the device immediately if there is any pain or discomfort.
- Keep fingers, hair, or other body parts away from the shaft and back of the vibrating head to avoid pinching.
- Do not place any object in the vent of the device.
- Do not immerse the device in water or allow water to enter the vent.
- Avoid dropping the device or using it improperly.
- Charge only with a DC 5V input (adapter not included).
- Inspect the device and battery carefully before each use.
- Do not modify or alter the device in any way.
- Never operate or charge the device unattended.

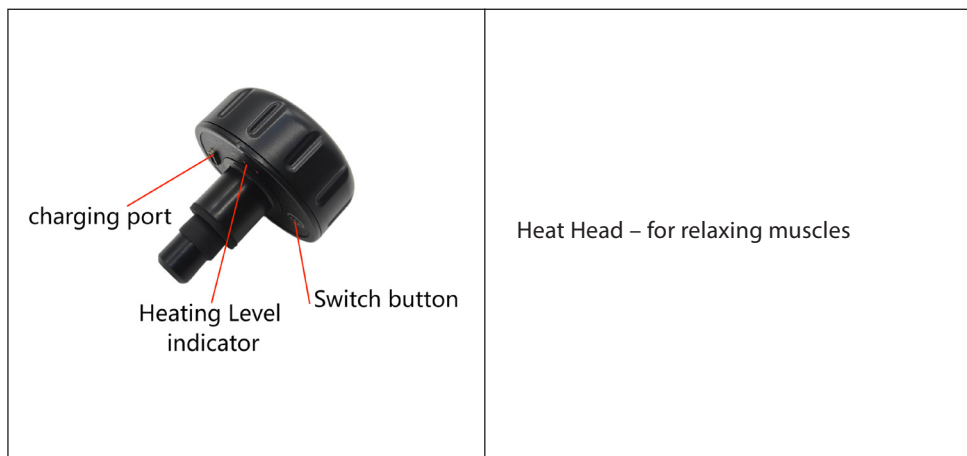
WARNINGS AND REMINDERS

- Do not immerse the device in water. Keep it away from liquids or heat sources.
- Keep the vent free from dust and debris.
- Do not remove the screws or attempt to disassemble the device.
- Do not use the device continuously for more than 30 minutes. Allow it to rest for 15 minutes before using it again.
- Remove the USB cable after charging or before use. (This product comes with a USB cable only; adapter not included.)
- For adults only. If you have any injuries, do not use this device. Please consult your doctor before use.

PART NAME & FUNCTIONS APPEARANCE CONSTRUCTION

 Massage Head Digital display screen Switch/ Mode Key Handle Charging port		 Level indicator Level adjustment button Mode adjustment button
		 Charging port

	U head/for vertebrae		Small flat head/for all parts of the body
	Ball head/for large muscle groups		Small round head/ for joints
		Cold Head / for sore muscles.	

Heat Head – for relaxing muscles

OPERATION INSTRUCTIONS:

1. Long press the power button for a few seconds to turn on the device. The default power level is 1.
2. Short press the power button to adjust the power level from 1 to 6. The levels cycle as follows:
1 → 2 → 3 → 4 → 5 → 6 → 1.
3. Press the mode button to select your preferred mode. The default is Mode 1.
 - P2: Mode 2
 - P3: Mode 3

NOTE:

- The device will automatically shut down after 15 minutes of use.
- The vibration head is removable and replaceable. To remove it, simply pull it out with a firm grip.
- Short press the power button to display the current battery level when the device is turned off and not charging.

USING THE HEAT HEAD

1. Press and hold the power button for 3 seconds to turn the device on. The red LED light ring will illuminate.
2. Press the power button once to cycle through the heating levels [Level 1 → Level 2 → Level 3 → Level 1].
3. To turn off the device, press and hold the power button for 3 seconds.
4. The device will automatically switch off after approximately 15 minutes.

MAINTENANCE, CLEANING, AND STORAGE:

- Clean the device with a slightly damp towel, then wipe it dry with a soft cloth.
- For storage, keep the device in a safe place or place it back in the supplied carton

when traveling.

CHARGING

1. Massager

- Before using the device for the first time, charge it fully by connecting it to a DC 5V input (adapter not included).
- When the device is charging, the LED screen will display the battery level to indicate the charging status.
- When the LED screen shows "99," the battery is fully charged.
- The average operating time is 75 minutes, depending on the speed setting and the pressure applied during use.

2. Heat Head

- Before first use, fully charge the heat head using a DC 5V input (adapter not included).
- The red light flashes while charging and remains solid when fully charged.
- Charging takes approximately 2–3 hours. After a full charge, the usage time is about 1.5–2 hours.
- For safety reasons, do not use the heat head while it is charging.
- To preserve battery life, avoid letting the battery fully discharge before recharging.

TECHNICAL DATA

Article No. 11088811

Input: DC 5V, 1A

20 Watt

1200mAh Li-ion battery +

300mAh Li-ion battery

PRODUCER

Epiq ApS

Krøyer Kielbergs Vej 3, 3. th

8660 Skanderborg, Denmark

www.epiq.dk

We cannot be held responsible for any printing errors.