salling

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Vi tager forbehold for trykfejl.

INTRODUCTION

Thank you for choosing a product from Salling. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept out of reach of children.
- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity,

physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.

PRODUCT DISPOSAL

Used product should not be disposed as household waste. According to legislation, the product should be disposed of properly so that the materials can be recycled in an environmentally friendly manner. Further information is available from the

responsible local authority or local disposal company. The product is equipped with the following logo:

PRODUCT FUNCTION

- A. The powerful pulsed water flow thoroughly cleans bacteria, tooth stains, tartar, and food debris from all areas of the oral cavity.
- B. The ultrafine water flow effectively reaches and cleans the gum line and deep spaces between teeth, promoting gum health.
- C. It efficiently removes stubborn dental plaque from tooth surfaces, offering a significantly better effect than traditional dental floss.
- D. Quick oral cleaning: Whether at home or dining out, it is the most effective way to rapidly restore oral cleanliness.
- E. Excellent for peripheral cleaning of orthodontic braces, dental implants, and veneers.
- F. Encourages healthy daily interdental cleaning habits, keeping teeth fresh and clean at all times.

POINTS FOR ATTENTION

DANGER

- A. Please read this manual carefully before use to avoid injury or damage caused by improper operation.
- B. Do not disassemble the product without authorization to prevent product damage or electric shock.
- C. Do not use water above 40°C to flush teeth in the water tank to avoid burns.
- D. If the product accidentally falls into water while charging, first unplug the power plug, then remove the product, and only recharge it after the water stains are completely dry.
- E. If the power cord or plug is damaged, please replace the power adapter.
- F. Avoid placing the product in high-temperature environments. Do not expose the product to fire, as the battery may explode.

WARNING

- A. Do not use the product for any purpose other than oral cleaning.
- B. Do not press the nozzle on the gums or teeth to avoid damaging the gums.
- C. Do not immerse the product in water to prevent product failure.
- D. Do not use the product in high-temperature and humid environments.
- E. If the power cord is damaged, please do not continue to use it.

ATTENTION

- A. Do not use this product without installing the nozzle.
- B. Some functional mouthwashes may be corrosive. Please check the composition of mouthwashes before use. Do not use them indiscriminately.
- C. If children, elderly individuals, or other operators wish to use this product, ensure that someone supervises and instructs the user.

- D. Keep the product out of the reach of children and infants. Do not use it as a toy.
- E. Patients with periodontal disease, those who have had oral surgery within the past 2 months, or individuals with oral issues should consult a doctor before use and only proceed after confirmation.
- F. Do not use the strong setting to wash damaged gums continuously, as it may cause bleeding or gum inflammation.
- G. Adjust the intensity of the pulsed water based on your oral condition in real-time.
- H. Do not immerse the entire unit in water, as this may cause product failure or complete damage.
- I. If the product is not used for an extended period, ensure the water tank is dry and free of water before storing to prevent bacterial growth.

KIND REMINDER:


When using this product, tilt it slightly.

Do not stand it upside down or tilt it too far to prevent liquid leakage.


USE OF RECHARGEABLE BATTERIES

- Before using this product for the first time, please fully charge it.
- During subsequent uses, if the three mode indicator lights flash simultaneously, it indicates that the product needs to be charged.
- When charging, use the built-in USB charging cable to connect the product to a power supply with an output voltage of 5V and a current of no more than 2A. While charging, the three mode indicator lights will flash. Once fully charged, the indicator lights will remain on.

STRUCTURE OF DENTAL PUNCH

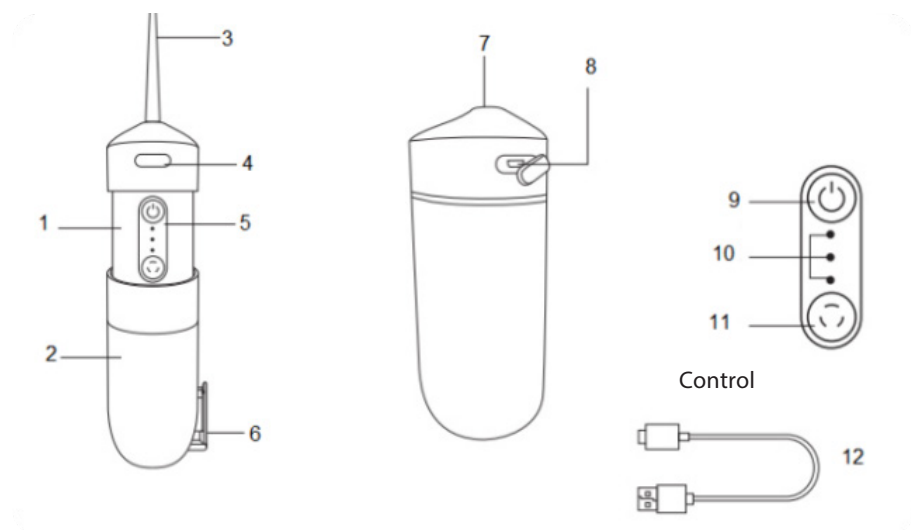1. Grip Handle
2. Water Tank
3. Nozzle
4. Nozzle Release Button
5. Control Panel
6. Water Supply Cover
7. Nozzle Mounting Hole
8. Charging Interface
9. Start/Close Gear Adjustment (Third Gear)
10. Gear Indicator Light
11. Mode Switch Key
12. Charging Cable

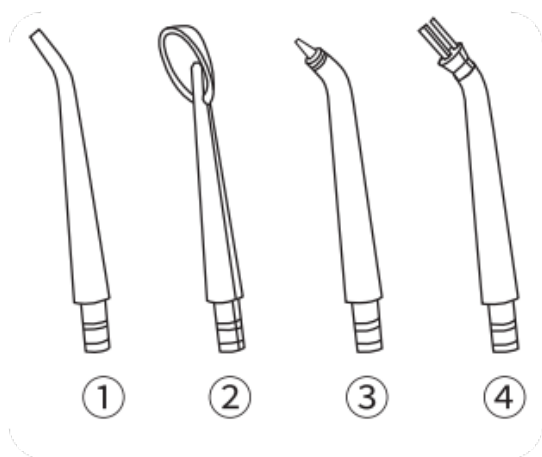1. Standard Nozzles: Suitable for daily cleaning.
2. Tongue Coating Nozzle: Specially designed for cleaning tongue coating.
3. Periodontal Pocket Nozzle: Specially designed for cleaning periodontal pockets.
4. Orthodontic Nozzle: Suitable for cleaning dental bands.

OPERATION METHOD:



BEFORE USING THIS PRODUCT:

- When using this product for oral cleaning for the first time, you may experience gum bleeding even if your gums are healthy. This is because your gums are being stimulated for the first time. Generally, bleeding will stop after 1-2 weeks of use.
- If gum bleeding persists for more than 2 weeks, it may indicate an abnormal condition (such as periodontal disease). In this case, please stop using the product and consult a dentist.

USAGE INSTRUCTIONS:

- Please install the nozzle before use.
- It is recommended to use running tap water or warm water. If the water temperature feels too low, you can use warm water that does not exceed 40°C.
- Do not turn on the switch before placing the nozzle in your mouth.

1. Open the Water Supply Cover:
 - Open the water supply cover to begin the setup.
2. Pull the Body to Its Intended Position:
 - Pull the body up to the correct position.
 - Note: When the water supply cover is closed, it may be difficult to pull out the body.
 - Pull the body out to the appropriate position and slide the limit buckle of the water tank into the lower limit slot of the body.
3. Insert the Nozzle:
 - Hold the body and insert the nozzle into the nozzle insertion hole.
 - Insert it firmly until you hear a click.
4. Inject Water:
 - Inject water into the water tank.
 - Do not inject chemicals such as salt, toothpaste, or mouthwash, as this could cause damage to the product.
 - To ensure continuous use, please turn off the switch and add water when needed.
5. Close the Water Supply Cover:
 - After adding water, close the water supply cover tightly to prevent leakage.
6. Place the Nozzle into Your Mouth:
 - Make sure the nozzle is positioned correctly in your mouth.
 - If the product is tilted excessively, the suction pipe inside the water tank may detach from the water surface, causing the water flow to stop.
7. Turn on the Switch and Select the Mode:
 - Short press the power button to start the device.
 - The default mode is "Soft Mode," and the soft mode indicator light will turn on.
 - Under full water conditions, the usage time is approximately 40 seconds.

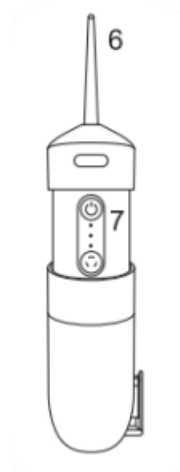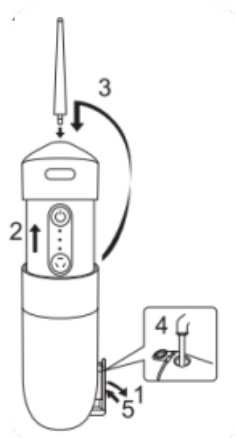SELECTION MODE (INTENSITY):

- Short press the mode key to switch between Standard Mode, Pulse Mode, and Soft Mode.
- Select your desired mode, and the corresponding mode indicator light will light up.

IMPORTANT NOTES:

- Do not press the tip of the nozzle too tightly against the teeth or gums.
- Only operate the product with water in the water tank (except for drainage purposes). Operating without water may cause a malfunction.



PROMPT:

- Hold the product vertically with the nozzle facing upwards.
- Lean your face over the sink.
- Place the nozzle into your mouth.
- Slightly close your mouth to prevent splashing.

AFTER USING THE PRODUCT:

After use, there is a risk of bacterial growth, so it should be thoroughly cleaned and dried before storage.

1. Open the water supply cover and pour out the water from the water tank.
2. Turn on the switch to drain the water inside the product.
- Run the product for 2 to 3 seconds until no more water is discharged.
3. Turn off the switch.
- Short press the on button to shut down, with automatic timed shutdown after 2 minutes.
4. Hold the body and pull out the nozzle.
5. Store the body.
- Press on the body until it clicks into place.
- Be careful not to pinch your fingers on the upper part of the water tank.
6. Close the water supply cover.

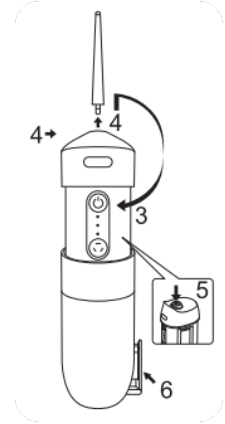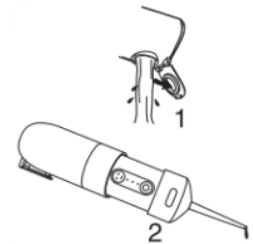TIPS:


- If the product needs to be stored in a bag or similar items after use, please wrap it in a clean towel or another absorbent material to prevent residual moisture from wetting the storage bag.
- To prevent bacterial growth, the product should be removed from the bag as soon as possible and allowed to dry completely.

CORRECT DENTAL CLEANING METHOD:

1. Choose the gear that suits you. It is recommended to start with the lowest gear and gradually adjust to the setting that works best for you.
2. When cleaning your teeth, insert the nozzle into your mouth, aligning it with the gaps or gums. Slightly open your mouth to ensure the water can flow out smoothly during the cleaning process.
3. After selecting your preferred gear, start from the innermost point (such as position 1 or 4) and rinse the gums in the same direction until you reach the gum line of the last

tooth.

4. After flipping the nozzle, align it with the gums and continue rinsing from the inside (such as positions 2 or 3) until you reach the gum line of the last tooth.
5. Repeat the same flushing sequence for the other row of teeth.

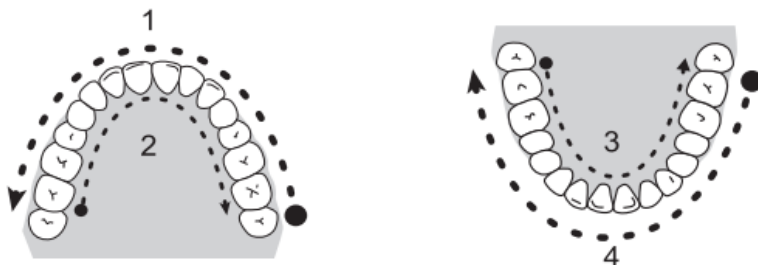WARM REMINDER FOR PRODUCT USE:


Dental experts recommend that, for first-time users, you should start with the gentle mode for a week. After getting used to it, you can switch to the normal mode by pressing the switch button. While using the product, please note the following:

- A. The direction of water flow should be perpendicular to the teeth or gums.
- B. Move slowly along the teeth, cleaning them in sequence.
- C. Align the nozzle perpendicular to the gum.
- D. Ensure thorough cleaning of teeth, braces, crowns, and dentures.

For first-time use, avoid spraying water directly into the periodontal pockets for extended periods to prevent gum discomfort.

√ Correct
demonstration

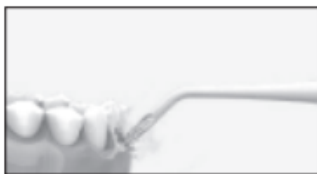× Error
demonstration

TIPS FOR PUNCHING TEETH:

- Always place the nozzle into your mouth before turning on the switch to prevent water from splashing everywhere.
- Gradually increase the gear from weak to strong during use. This allows your gums to adapt to the water pressure and helps you find a suitable gear that is comfortable for your gums.

ROUTINE MAINTENANCE:

- A. If you need to clean this product, rinse it with water or wipe it clean with a damp cloth. Then, use a dry cloth or other dry tissue to wipe off any remaining water stains on the surface.
- B. It is strictly prohibited to immerse the product body in water.
- C. Do not use hot water exceeding 50°C to clean this product.
- D. Avoid bending, stretching, or twisting the nozzle water pipe.
- E. It is recommended to replace the nozzle every 3 months.
- F. If you do not use this product for an extended period, empty the water tank and store the product in a cool, dry place.

TROUBLESHOOTING

Common Problems	Possible Causes	Solutions
Unable to work or power on	Newly purchased or not used for more than three months	Please charge for at least 4 hours before use
The water pressure is very low during operation	The nozzle is clogged or deformed	Replace the nozzle
It can only be used for a short time while charging	Insufficient charging time	Please Charge for more than 4 hours before use
When charging, the charging indicator does not light up or flashes	The power supply is not connected properly	Please ensure the power is connected correctly
	Adapter (charging head) mismatch	Replace the adapter
Do not spray water when in use	Newly purchased or unused for a period of time	1. Pull out the nozzle from the head 2. Blow hard at the nozzle to remove any blockages

TECHNICAL DATA

Article No. 10958380
5V/1A, 5 Watt
Capacity: 200 ml
Charging time: 3.5 hours
Usage time: 40 minutes
800mAh Li-ion battery

PRODUCER

Epiq ApS
Krøyer Kielbergs Vej 3, 3. th
8660 Skanderborg, Denmark
www.epiq.dk

We cannot be held responsible for any printing errors.