&me.

INSTRUCTION MANUAL MASSAGE GUN

10958218

INTRODUKTION

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Vi tager forbehold for trykfejl.

INTRODUCTION

Thank you for choosing a product from &me. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- The product engine/electronics, cord or plug cannot be submerged in water, and no water may get in contact with the engine/electronics.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept out of reach of children.

- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity, physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.

PRODUCT DISPOSAL

Used product should not be disposed as household waste. According to legislation, the product should be disposed of properly so that the materials can be recycled in an environmentally

friendly manner. Further information is available from the responsible local authority or local disposal company. The product is equipped with the following logo:


Warning: To minimize the risk of electric shock, fire, and personal injury, ensure that this product is used in strict accordance with the following instructions.

- For adults only.
- Use only on dry, clean skin surfaces of the body; do not use through clothing. Gently press and move it on the skin. Duration on each part should be approximately 60 seconds.
- Use this device only on soft tissues of the body without causing pain or discomfort. Avoid use on the head or any hard or bony part of the body.
- Only use the vibrating head that produces the best results.
- Bruising may occur under moderate or normal pressure. Always check the area receiving physiotherapy.
- Stop use immediately if there is any sign of pain or discomfort.
- Keep fingers, hair, or other body parts away from the shaft and back of the vibrating head to prevent pinching.
- Do not place any object in the vent of the device.
- Do not immerse it in water or allow water to enter the device's vent.
- Avoid dropping the device or using it improperly.
- Charge the device with DC5V 1A input only.
- Carefully inspect the device and battery before each use.
- Do not change or modify the device in any way.
- Never operate or charge the device in an unattended manner.

WARNINGS AND REMINDERS

- Avoid immersing the device in water. Keep it away from liquids and heat sources. Ensure the vent remains free from dust and debris.
- Do not attempt to remove screws or disassemble the device.
- Limit usage to no more than half an hour continuously. Pause the device for 15 minutes before resuming.
- Unplug the charger (this product comes with a USB cable only; charger not included) from the massager after charging or before use.
- This device is intended for adults only. If you have injuries, refrain from using this device, and consult with your doctor before use..

KNOW YOUR MESSAGE GUN

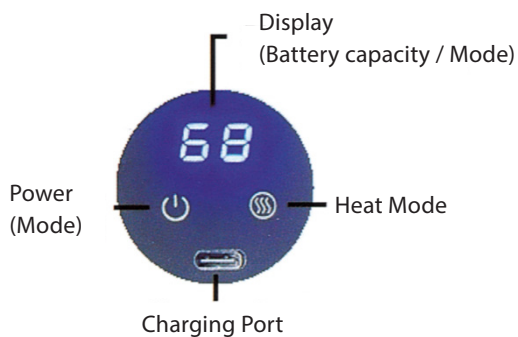Select Speed Display

Heat Mode is not activated

Heat Mode is activated

Speed at Level 1: 2800-2900rpm
 Speed at Level 2: 2900-3000rpm
 Speed at Level 3: 3000-3200rpm
 Speed at Level 4: 3200-3400rpm
 Speed at Level 5: 3400-3600rpm
 Speed at Level 6: 3600-3800rpm

Cold Head / for muscle inflammation.

Heat Head / for soothes aching muscle

U head/for vertebrae

Small round head/ for joints

Ball head/for large muscle groups

Small flat head/for all parts of the body

HOW TO USE

- Choose a suitable massage head and plug it into the massage gun.
- Long-press the Power button for 1-2 seconds to turn it on.
- Short press to switch between 6 different levels of vibration.
- Auto Shut Off is set at around 15 minutes.

Notice: When the massage movement begins to slow down, stop using the massager, and recharge it using the supplied charger.

HOW TO USE THE HOT-HEAD

- Plug the Hot Head into your massage gun.
- Long-press the Power button for 1-2 seconds to turn on the appliance; the red light will illuminate on the Heat Mode button. Wait for 3 to 4 minutes for the Hot Head to heat up.
- Set the massage gun to the desired speed and place it on the desired area.

HOW TO USE THE COLD-HEAD

- Before using the Cold Head, place it in a freezer for about 2 hours for best results.
- Plug the Cold Head into your massager.
- Long-press the Power button for 1-2 seconds to turn on the appliance.
- Set the massage gun to the desired speed and place it on the desired area.

MAINTENANCE, CLEANING, AND STORAGE

Clean the appliance with a slightly damp towel, and then wipe it dry with a soft cloth for storage or place it in the supplied carton when you go traveling.

CHARGING

- Before using it for the first time, fully charge the appliance (approximately 3-4 hours) using a 5V 1A adaptor (not included).

- The battery capacitor indicator will turn on and rise as the battery is being charged.
- The battery capacitor indicator will reach “99” when the battery is fully recharged.

TECHNICAL DATA

Article No. 10958218

Input: DC 5V, 1A

25-45 Watt

2 x 1800 mAh Li Battery

PRODUCER

Epiq ApS

Krøyer Kielbergs Vej 3, 3. th

8660 Skanderborg, Denmark

www.epiq.dk

We cannot be held responsible for any printing errors.