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11298757

INTRODUKTION

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Vi tager forbehold for trykfejl.

INTRODUCTION

Thank you for choosing a product from Essential. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept out of reach of children.
- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity,

physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.



PRODUCT DISPOSAL

Used product should not be disposed as household waste. According to legislation, the product should be disposed of properly so that the materials can be recycled in an environmentally friendly manner. Further information is available from the responsible local authority or local disposal company. The product

is equipped with the following logo:


PRODUCT DESCRIPTION

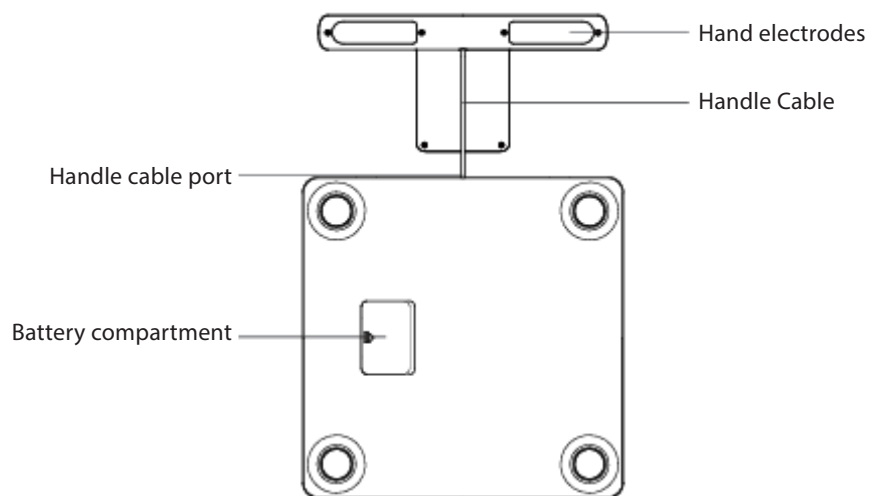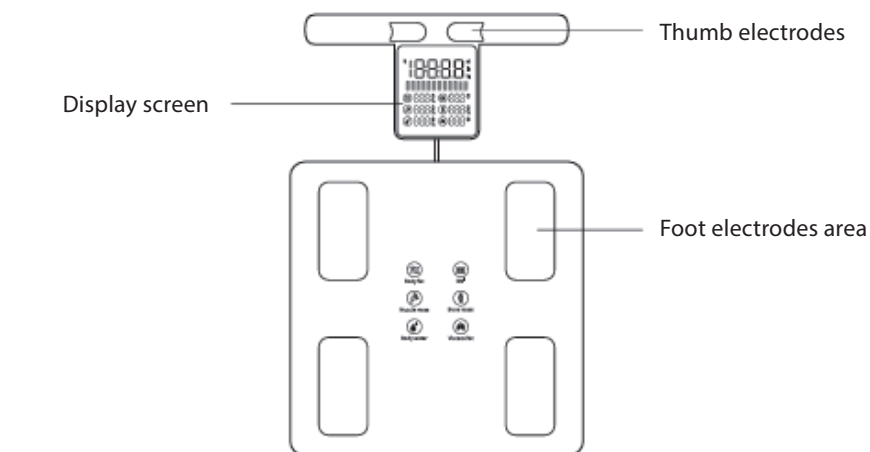Product Back Demonstration

INSTRUCTIONS FOR USE

Screen Display

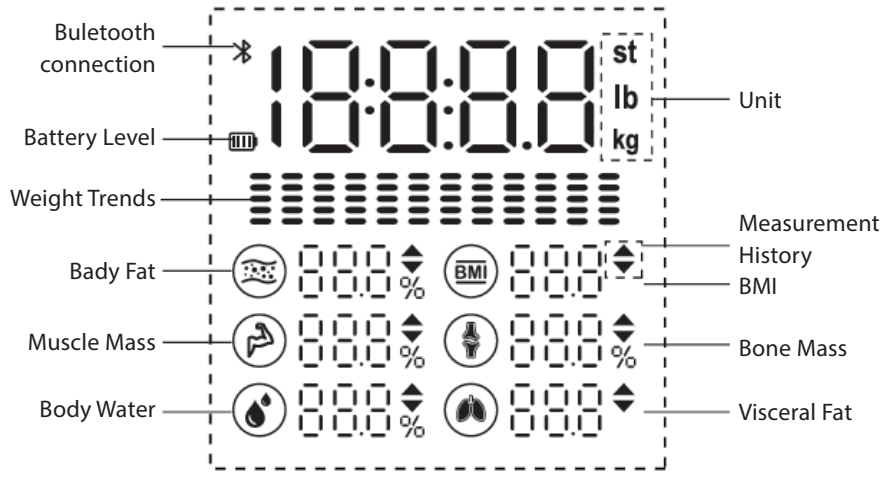DISPLAY SCREEN ICON GUIDE

	Body fat : Illuminates when displaying body fat percentage measurement results.
	Muscle icon: Illuminates when displaying muscle mass measurement results.
	Body water icon: Illuminates when displaying body water measurement results.
	BMI icon: Illuminates when displaying BMI measurement results.
	Bone mass icon: Illuminates when displaying bone mass measurement results.
	Visceral fat level icon: Illuminates when displaying visceral fat level measurement results.

 Normal working	No display Change battery	 No wireless connection	---- Return handle
 Change battery	 Overload	 Power cycle	 Calibrating

INITIAL SETUP

Step 1: Download app and create an account

From Apple App Store (ios) or Google Play (Android), search for “Feelfit” and download the app to your smart phone. You can also scan the QR code to download.

Scan QR code
and download “Feelfit” app


Register an account


App icon

Step 2: Installing batteries (3xAAA not included)

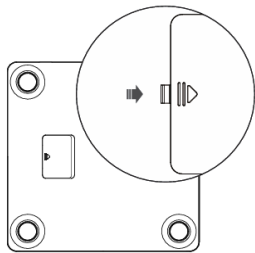1. Slide the lid in the direction indicated by the arrow.

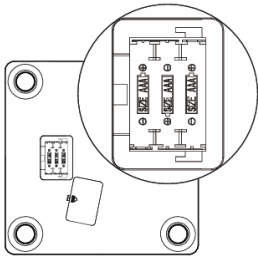2. Insert the battery and then close the cover.

Step 3: Pair scale with smartphone.

Turn on the Bluetooth on your phone. The product will start to automatically pair. At the same time, the name of the device will be displayed.

Note: For Android 6.0 or above, please also turn on Location.


MEASUREMENT STEPS

Preparatory work

1. Place the scale on a clean, flat, hard surface. Wait until the display reads 0, then connect the APP during this period.
2. Connect App to Create User Profile.
3. Try wiping them lightly with a soft towel that has a little water on it. This will help get a better reading.

WEIGHT MEASUREMENT

Stand barefoot on the scale with your body upright. Grip the handles firmly with palms facing upward. Once the weight reading stabilizes, the flashing weight indicator signals completion.

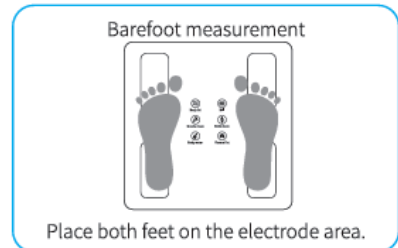BODY FAT SCALE MEASUREMENT

- Hand Demonstration Image

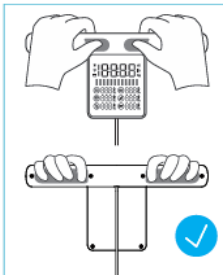
Fig. 1

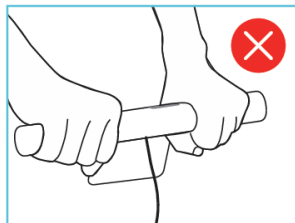
Fig. 2



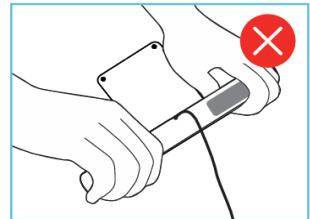
Fig. 3


1. Place your thumb on the oval electrode and your four fingers around the elongated electrode, with the handle connector wire pointing downward.
 2. Hands should not be clenched.
 3. The handle connector wire should not face upwards.
- Body Demonstration Image

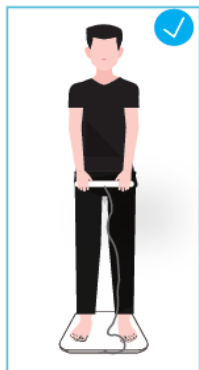
1. Stand upright with arms naturally hanging down at an angle from the body, ensuring arms do not touch the waist. Undgå at bøje hænderne udad.

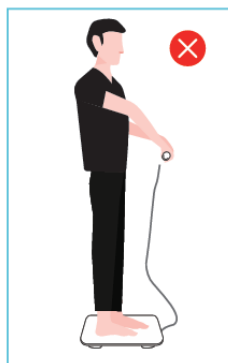
2. Avoid bending the hands outward.

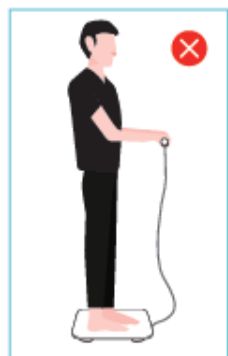3. Do not extend the hands horizontally for measurement.

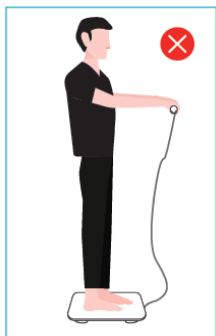4. Avoid bending the hands during measurement.

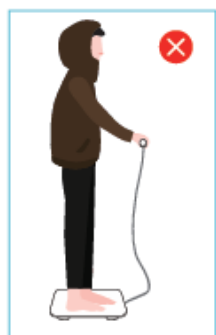5. Wear formfitting clothing only during the measurement process.

• Foot Placement Illustration

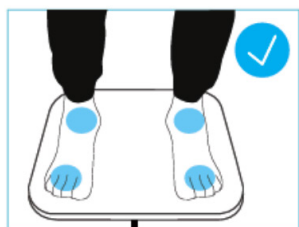
Fig. 1

Fig. 2

Fig. 3

1. Stand with feet parallel on the scale electrodes, ensuring inner thighs do not touch each other.
2. The heels of both feet should not touch each other.
3. The feet should not be placed in a heel-to-heel stance.

CAUTION

1. Do not use if you have a pacemaker, or any other internal medical devices in your body.

2. Do not use if you are pregnant.
3. Do not use when your body is wet.
4. Do not use rechargeable 1.2V AAA batteries.
5. Do not use in damp environments such as bathrooms; the scale is not waterproof.

TROUBLESHOOTING

Problem	Cause	Solutions
	Low power	Please replace and use new 1.5V AAA battery.
	Overload	The scale within 400lbs/180kg.
	Weight detected after stepping off	When this icon scrolls, place the handle back on the scale and position it on a clean, flat and firm surface. Wait for the scale to complete calibration.
No Display	Auto-off, or no power	Auto-off: Step on then off to turn on the scale; No power: Please change and use new 1.5V AAA battery.
Display Weight Only	Unable to determine body composition	User Setup: Please connect the app and create a user profile before use Age Restriction: Ensure the age setting is 10 years or older. Posture Requirements: Measure according to the standard measurement posture.
Unable to Sync with App	Lost connection, or scale turns off	Reconnect the scale: For iPhone, authorize Bluetooth; for Android, enable location. Click the icon on the top-right corner and tap "My Device" > "+" Add device" on the app.
Weighing Deviation	Uneven floor surface, contaminated scale feet, undetected scale body weight	Unstable surface: Place the measuring device on a flat, solid surface. Dirty base: Remove debris from the device surface and base. Reset device: Remove the battery, wait 3 seconds, then reinsert it. Place the device on the ground and wait for the display to turn off before normal use. After measurement, return the handle to the mobile unit to view the full report.

TECHNICAL DATA

Article No. 11298757

Battery: 3xAAA (not included)

PRODUCER

Epiq ApS

Krøyer Kielbergs Vej 3, 3. th

8660 Skanderborg, Denmark

www.epiq.dk

We cannot be held responsible for any printing errors.