BRUGSANVISNING IPL HÅRFJERNER

10958382

INTRODUKTION

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Vi tager forbehold for trykfejl.

INTRODUCTION

Thank you for choosing a product from Essential. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- The product may only be connected to 100-240V~, 50/60Hz.
- The product engine/electronics, cord or plug cannot be submerged in water, and no water may get in contact with the engine/electronics.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept

out of reach of children.

- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity, physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.

PRODUCT DISPOSAL

Used product should not be disposed as household waste. According to legislation, the product should be disposed of


properly so that the materials can be recycled in an environmentally friendly manner. Further information is available from the responsible local authority or local disposal company. The product is equipped with the following logo:

WARNING


- In order to achieve the best results and avoid malfunctions, please ensure you carefully read the entire user manual before using the device.
-  Always wear goggles when using the device. Do not use it directly on the eyes, as the strong light emitted by this device can cause eye damage.
-  Never use the appliance near water (e.g., over a filled wash basin, bathtub, or in the shower).
- Do not use near children, as the device emits strong light that can easily reach the eyes.
- Regularly check and clean the irradiation port to avoid burns and malfunctions during use.
- Do not use the device when the ventilation port is blocked to prevent smoke or fire.
- Do not direct the irradiation port towards black objects, as this can prevent the machine from functioning properly.
- Do not expose clothing, hair, or other items to the irradiation port to prevent burning, discoloration, or damage.
- Do not use on pets to avoid injury.
- If the skin experiences slight warmth and redness after use, this is a normal reaction. You can use the cooling compress function to alleviate any discomfort.
- When removing hair from the back, it is recommended to have someone assist you to avoid risks.
- After turning on the product, never leave it unattended. Turn it off promptly after use.
- If there is a temperature difference of more than 20°C between the operating environment (from a low-temperature to a high-temperature environment or vice versa), turn off the device and allow it to stand for 2 hours before use.
- Avoid direct ultraviolet radiation on the irradiated area within 24 hours after treatment.
- If the skin of the treated area is not dry and clean, it will affect the effectiveness of the pulsed light and may cause the product not to flash or accidentally flash.
- If the built-in fan does not run (no buzzing sound is heard) after turning on the power, do not use the device to prevent burning, smoking, or fire.
- Do not use the IPL Handset if you are under the age of 16.
- Do not use the IPL Handset if you are menstruating, pregnant, or lactating.
- Do not use the IPL Handset on tattoos.
- Do not use the IPL Handset on broken or burned skin.
- Do not use the IPL Handset on the same area more than 3 times in a single session.

PARTS DESCRIPTION

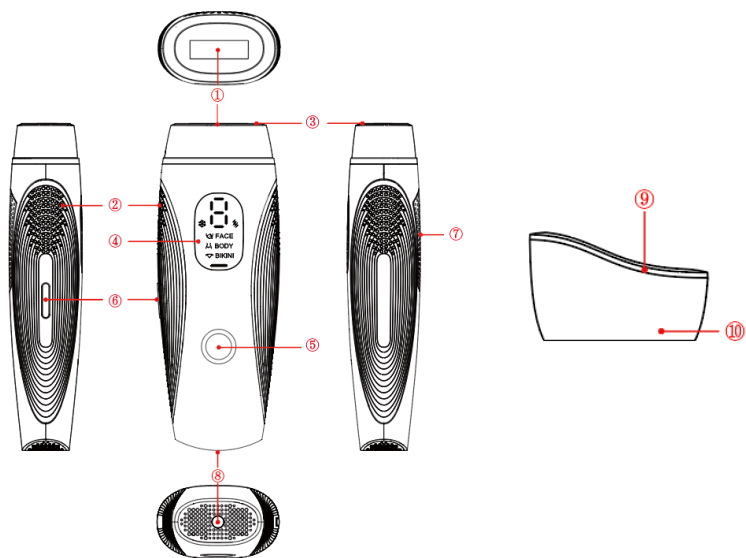- | | |
|------------------------|-----------------------|
| 1. Light Outlet | 6. Power/Gear Button |
| 2. Air Outlet | 7. Air Inlet |
| 3. Cold Compress Panel | 8. DC Socket |
| 4. LED Display | 9. Disinfection Lamp |
| 5. Flash Button | 10. Disinfection Base |

LED DISPLAY AND MODE INTRODUCTION

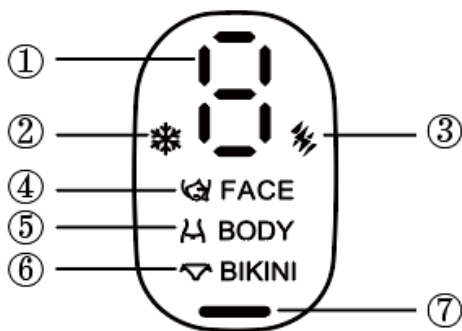
- Levels Icon (1-9):
 - This IPL hair removal device has 9 intensity levels.
 - Short press the Power Button/Gear Button to increase or decrease levels cyclically (1 → 9 → 1)


2. Cooling System Icon  : Always lights up after startup.
3. Flashing Mode Icon  :
 - Press and hold the Flash Button for 3 seconds to switch between Manual Flashing Mode and Automatic Flashing Mode.
 - Repeat to toggle back to Manual Flashing Mode.
4. Face light indicator icon  **FACE** :
 - Always lights up when the product setting level is 1 to 6.
 - When the level sets exceed 6, the indicator light turns off.
5. Body light indicator icon  **BODY** :
 - Always lights up when the product setting level is 1 to 9.
6. Bikini light indicator icon  **BIKINI** :
 - Always lights up when the product setting level is 1 to 3.
 - When the level sets exceed 3, the indicator light off.
7. Induction indicator light: After the machine is powered on, the indicator light enters a flashing state. When the light outlet is close to the surface of the skin, the indicator light remains on.

BEFORE USING THE APPLIANCE

1. Thoroughly clean your skin, ensuring it is free of hair, completely dry, and free of oily residues such as perfume, lotion, deodorant, or sunscreen.
2. Ensure that the irradiation port is free of dust, hair, or any dirt, and confirm there are no damages or abnormalities. Clean it thoroughly if necessary to prevent product malfunction, skin irritation, or other issues.
3. Check the operating environment temperature of the device. It should only be used in a room with a temperature between 5°C and 40°C.
4. Adjust the irradiation intensity based on your skin condition. For first-time use, always start at level 1.
5. Do not use this device on difficult or sensitive areas, such as ankles or bones.
6. Avoid using the device on areas with burn marks, dark pigmentation, moles, frostbite, sunburn, eyebrows, or around the eyes.

OPERATION INSTRUCTIONS

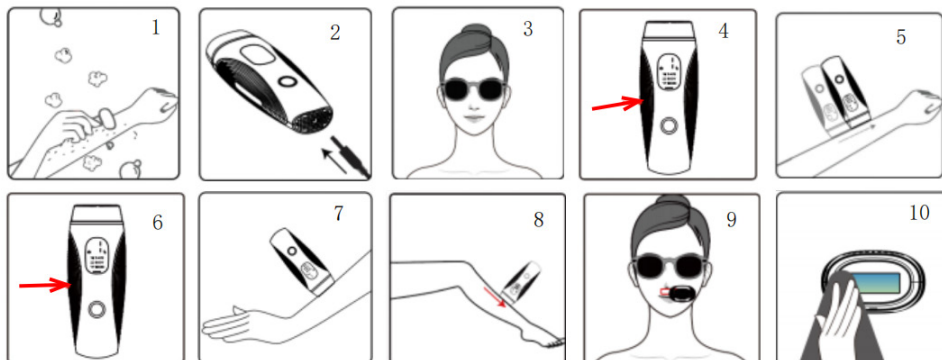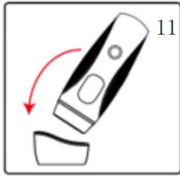
1. Remove hair from the target area, clean the skin, and ensure it is completely dry before use. (See Figure 1).
2. Connect the IPL hair removal device to the power adapter. (See Figure 2).
3. Put on the protective goggles included with the product. (See Figure 3).
4. Press and hold the Power Button/Gear Button to turn on the device. You will hear a sound, and the relevant parameters will display on the screen. The default gear level is set to 1. (See Figure 4).
5. Align the device's Light Outlet with the desired area for treatment. When the induction prompt light turns on, press the Flash Button to start. (See Figure 5).
 - To switch between Manual and Automatic Flashing Modes, press and hold the Flash Button for 3 seconds.
 - For areas near the eyes, it is recommended to use Manual Mode for better control.
6. Use the Gear Button to adjust the intensity level. Short press to gradually increase the gear level (1 is the lowest, 9 is the highest). (See Figure 6).
7. Ensure the Light Outlet is in close vertical contact with the skin. When the induction prompt light is on, press the Flash Button to emit a light pulse. (See Figure 7).
 - Note: The Flash Button is used to activate flash pulses and supports both Manual and Automatic Flashing Modes.
8. Avoid flashing repeatedly on the same area, as this can cause skin damage or burns. Slide the device across the treatment area while flashing. (See Figure 8). (See Figure 9).
9. Once finished, press and hold the Power Button for 2 seconds to turn off the device. Disconnect the power adapter and allow the device to cool completely before cleaning. (See Figure 10).
10. Place the appliance on the base for disinfection. You will see the blue light turn on, indicating the appliance is disinfecting. (See Figure 11).
 - The disinfection base must be plugged into the USB Type-C charging wire. Inspect the lamp head of the hair removal device and ensure it is free of dirt.
 - Ensure the disinfection time is more than 30 minutes to achieve the sterilization effect.

IMPORTANT NOTES

1. Do not use this device on naturally black skin.
2. Avoid direct sunlight exposure before performing hair removal.
3. Shave the hair before using this device. (It is recommended to avoid chemical hair removal methods).
4. Start at Level 1 when using the device for the first time.
5. Always wear the provided goggles while using the device.

6. Do not use this device in the dark.
7. Avoid using the device under direct sunlight to prevent damage or malfunction.
8. Handle the device carefully; do not strike or shake it violently.

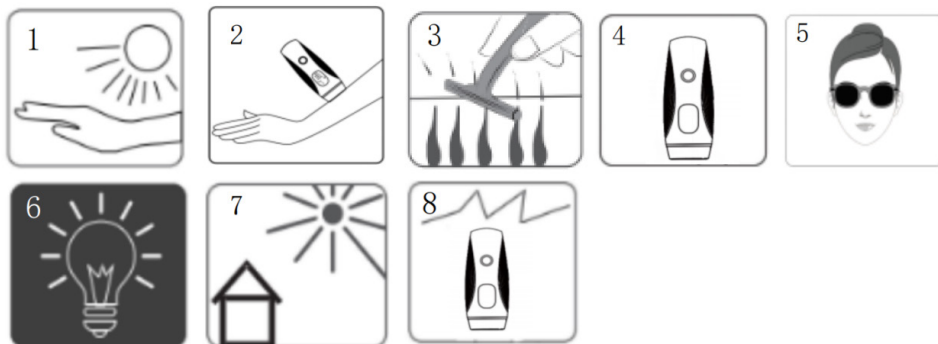USAGE INTERVAL

- Weeks 1-4: Use 1-2 times per week.
- Weeks 5-8: Use once a week.
- After Week 8: Use once a month or as needed.

For optimal results, use the device based on hair growth patterns. For individuals with particularly fast hair growth, the intensity (gear level) and frequency of use can be increased accordingly. The ideal usage may vary from person to person.

CLEANING AND MAINTENANCE

- To ensure optimal performance and extend the product's lifespan, clean the device before and after each use to prevent the accumulation of hair and dust. Improper cleaning may affect the device's performance.
- Turn off the device and unplug the power adapter. Allow it to cool completely before cleaning.
- Do not clean the device or its accessories under running water.
- Avoid using steel wool, abrasive cleaners, or corrosive liquids when cleaning the product.
- Note: If the transparent window becomes difficult to clean, discontinue use of the device.
- If necessary, use a clean, soft cloth to wipe the exterior of the product.

STORAGE

- Inspect the product, particularly the transparent window and body components, to ensure they are clean and free of debris.
- When not in use, store the appliance in a dry, well-ventilated location, out of reach of children.

- Place the power adapter near the appliance, avoiding wrapping the cord around it.
- Ensure the hair removal device is stored safely to maintain its longevity and effectiveness.

TECHNICAL DATA

Article No. 10958382
Input: DC12V, 3A
100V-240V, 50/60Hz
36 Watt

PRODUCER

Epiq ApS
Krøyer Kielbergs Vej 3, 3. th
8660 Skanderborg, Denmark
www.epiq.dk

We cannot be held responsible for any printing errors.