Overschiestraat 63
1062 XD Amsterdam
The Netherlands
info@optismile.eu

www.optismile.eu

Art. No.: 3031 / Batch no. 22626

Waterflosser Manual

NL • Handleiding FR • Manuel DE • Bedienungsanleitung
SE • Manual DK • Brugsanvisning NO • Bruksanvisning
FI • Käyttöohje ES • Manual PT • Manual

EN


Thank you for purchasing our water flosser. This manual will provide instructions for effective and safe use of the product. Please read it carefully before use and store the manual in a safe place for future reference. If you ever pass the product on to others, this manual should be included.

1. Package content

- Water flosser
- 4 Nozzles
- USB-C cable (1 m)
- Manual

2. Technical data

Product model:	HX-902
Material:	ABS
Net weight:	249 g
Water tank capacity:	220 ml
Charging voltage:	5 V
Charging current:	1 A
Rated voltage:	3.7 V
Rated power:	4 W
Water pressure:	90-120 psi, 3 settings
Charging time:	2-6 hours
USB charging cable:	type C, 1 m
Waterproof:	IPX7

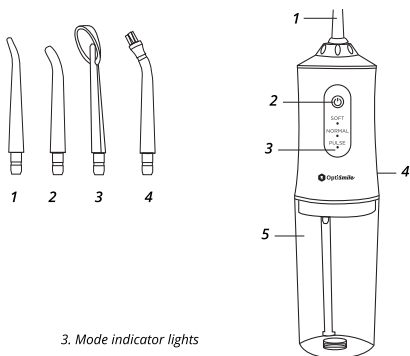3. Description of parts

- | | | |
|------------------|-----------------------|--------------------------|
| 1. Nozzle | 2. Power button | 3. Mode indicator lights |
| 4. Charging port | 5. Water tank (220ml) | |

Note: We always test the water jet of each product. That is why there may be some distilled water left inside. This is completely safe and will not effect your health.

1. **Standard nozzle:** for daily cleaning.

2. **Periodontal pocket nozzle:** to remove plaque and clean under the gumline. It features a soft rubber tip designed for low-pressure delivery of therapeutic rinses.

3. **Tongue scraping nozzle:** for gentle scraping and flushing. Removes bacteria and sulphur compounds that reside on the tongue and cause bad breath.

4. **Orthodontic nozzle:** tapered brush that helps remove plaque from braces and dentures.

4. Safety instructions

- For intended use only, as described in this manual.
- The product is safe to use by children from the age of 8 and by persons with reduced physical or mental capabilities, or with lack of experience or knowledge, but only under supervision and with taking into account any hazards.
- Maintenance should be done by or under the supervision of an adult.
- Before use, inspect the device and its parts. Do not use in case of damage.
- Do not let moisture come into contact with any electrical parts.
- Do not operate the device without a nozzle.
- Do not hold the nozzle too close to your teeth and gums.
- Remove any oral jewellery prior to use of this product.
- Do not use the device if you have sores, blisters or cuts in your mouth.
- When undergoing treatment, consult your dental professional prior to use.
- Stop using the product and contact your dentist, if you experience any discomfort or pain.
- This product is a personal care device and is not intended for use on multiple persons.
- Attachments not recommended by the manufacturer may cause a hazard.
- Do not attempt to repair, disassemble or alter the product yourself.

5. Operation - Charging

1. Remove the charging port cover.
2. Connect the water flosser to the mains using the provided USB-C cable and an adapter (not included).
3. The mode indicator lights start flashing alternately while charging and all three light up when the water flosser is fully charged.
4. Make sure the charging port cover is securely replaced before use.

Note: It takes 2-6 hours to fully charge the battery, but the water flosser can be used before that. When the mode indicator flashes during use, the water flosser needs to be charged.

Using the water flosser

1. Make sure the water flosser is fully charged (approx. 2-6 hours) before first use.
2. Place the nozzle of your choice onto the water flosser.
3. Remove the water tank and fill it with water (max 220 ml).
4. Short press the power button to turn on the water flosser. Long press to select one of these modes:
 - Soft: for sensitive teeth and children.
 - Normal: for deep removal of bacteria under the gums and between the teeth.
 - Pulse: for gum massage and cleaning.
5. After each use, rinse the parts carefully to remove any residue. Empty the water tank and wipe everything dry with a clean and soft cloth.

Note:

- Flossing is an addition to daily brushing and a part of a healthy dental routine.
- Make sure to clean all four quadrants of your mouth equally.

Smart timer

The water flosser will automatically shut off at the end of the dentist-recommended brushing time of two minutes.

6. Maintenance

- Before long-term storage, remove any remaining water to prevent bacteria from spreading.
- If necessary, clean with water or a mild non-abrasive detergent. Do not use water warmer than 40°C. Dry everything thoroughly with a clean and soft cloth.
- Store the product in a cool and dry place away from direct sunlight.

7. Troubleshooting

Problem	Cause of problem	Solution
The water flosser doesn't turn on.	The product is new or you haven't used it in 3 months.	Charge the water flosser for 24 hours.
The water flosser only works for a couple of minutes.	The battery died.	Please contact our customer service.
	The water flosser isn't charged.	Charge the water flosser for 24 hours.
The water pressure is very low.	The nozzle is congested.	Clean the nozzle or contact our customer service.
	The nozzle is bent.	
	There's dirt in the filter.	Clean the filter.
There isn't any water coming out of the water flosser.	The water tank is empty.	Fill the water tank with water.
	The water flosser is tilted too much.	Keep the water flosser straight.
Water drips from the product and water tank.	The water flosser is tilted too much, causing water to drain out of the vent or the water tank.	Keep the water flosser straight.

8. Disposal

Dispose of the article and the packaging materials in accordance with current local regulations. This product is subject to the provisions of the European Directive 2012/19/EU.

The crossed out wheeled bin symbol means that this product shall not be disposed of with normal household waste. Electronic and Electrical Equipment not included in the selective sorting process are potentially dangerous for the environment and human health due to the presence of hazardous substances. Please dispose of responsibly at an approved waste or recycling facility.

The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

Recycling helps to reduce the consumption of resources and reduce the environmental impact. If you wish to discard the packaging, comply to the environmental legislation in your country.

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