MASSAGE PRO

HEATING FOOT MASSAGE

BRUGERMANUAL

TECHNICAL PARAMETERS


Input voltage: AC220V 50Hz

Timer: 15 min


Rated power: 40W

ENG

MESSAGE DIAGRAM

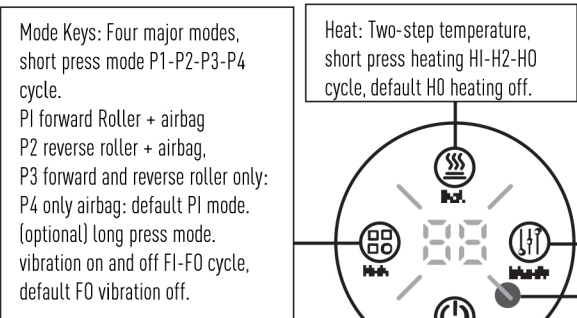*Notes: Whenever the operation button is pressed, the window indicator light will be on to indicate successful operation.

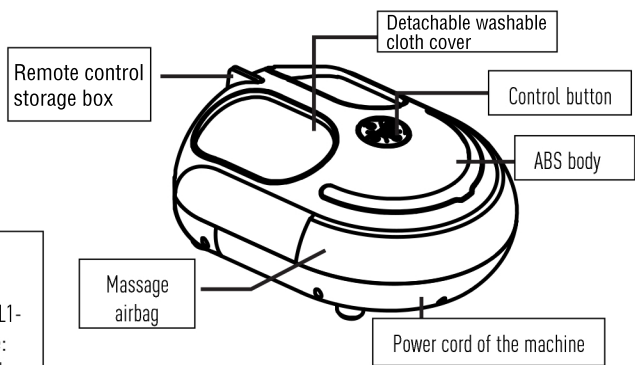APPLICABLE STANDARD

1. Pull the power cord plug out of the bottom slot and plug it into a power socket.
2. After connection to the power supply, the indicator light of the switch button on the key board is flashing
3. Press the switch button on the key board to start the machine.

SAFETY (CAUTIONS)

Please use the machine for not more than 30 minutes, otherwise it may cause pain to muscle and nerve.

1. The following users should consult a physician before using the product.
 - A user who is planted with a pacemaker or similar implantable medical device;
 - A user who suffers from heart disease;
 - A user who suffers from strain disease;
 - A user with fragile or injured bones;
 - Especially a user with injured back bone;
 - A user who just undergoes surgery;
 - A user who is pregnant or just gives birth;
 - A user who suffers from osteoporosis;
 - A user who suffers from blood circulation problems or thrombosis;
 - A user who suffers from skin infection;
 - A user who receives treatment or adjuvant treatment;
 - A user with abnormalities in the limbs;
 - A user who suffers spondylosis and vertebral injury caused by accident or disease.
2. If the cover is damaged or if there is mechanical and electronic failure, please do not use this product but contact the local seller or manufacturer.
3. Do not use this product in extremely humid and high-temperature rooms or use it beside heating objects (such as electric heating, stoves, etc.).
4. Before or during the use of this product, please do not apply medical or cosmetic creams or lotions to the product or your back.
5. The product can only be sent to the company's authorized repair center for testing. User shall not repair the product, otherwise the warranty will be invalidated
6. Do not vibrate or drop the product.
7. Do not use this product one hour before or after a meal.
8. Do not wind the power cord around the product or bend it between the upper back pads, otherwise it may be damaged.
9. Do not stand on the product or put a heavy object on it.
10. If the product fails, please unplug the power cord and immediately consult the manufacturer.
11. Electrical appliances are prone to causing risks after contact with water. Please do not use this product in a humid environment like bathroom or places near a faucet or swimming pool.
12. If the power cord is damaged, it must be replaced by the manufacturer or its service center or a person with relevant qualifications, so as to prevent risks.

SAFETY (WARNING)

Do not use power apply without specified voltage - The use of power supply with other voltage may conte accidents such as electric shock, fire, etc.
 Power outage - Please turn off the power switch and unplug power plug scared the service
 Machine failure - Please unplug the power plug immediately and contact the seller
 Machine failure - Please turn off the power and unplug the power plug
 Damage of power cord or plug - In order to prevent danger, replacement must be done by the manufacturer, or its maintenance department, or similar full-time personnel
 In operation - In order to prevent danger, it is forbidden to put hands into the cloth cover of the massager

MAINTENANCE METHOD

If there are stains on this machine after use, please wipe them with a piece of cloth wet with a neutral detergent. (Do not use chemicals such as petroleum essence, alcohol or diluent on the massager). Then, leave the machine for natural drying. (Please do not dry the machine in a intensified manner.)
 Warning: Before maintenance and cleaning, please turn off the power switch and unplug the power plug.
 Tips: This machine can be used continuously for about 30 minutes each time to achieve the ideal massage effect. Storage temperature: 0~40°C.

TROUBLESHOOTING

This product is controlled by a microcomputer, and is more sensitive to voltage fluctuations. If the voltage fluctuation exceeds 10% or -15% of the limited voltage, the following abnormalities may be caused.

POSSIBLE ABNORMALITIES & TROUBLESHOOTING

The machine it not working after the connection to the power supply.

1. Check whether the power cord of the machine is complete
2. Check whether the power plug is plugged or connected to power supply

The machine is not working after the connection to the power supply.
 If the temperature of the product and the engine is excessively high after long time operation, temperature protection will be automatically initiated Please turn off the power and leave the engine to dissipate heat for more than 30 to 50 minutes.
 The manual controller fails and the machine is not working during use.
 After the automatic timing is initiated, if the time is up, please press the button again to turn off the power and re-start the machine.
 If any abnormalities still exit after the troubleshooting, please contact the manufacturer or local seller.

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Stock Group A/S