MASSAGE PRO

Nakke og skuldermassage

BRUGERMANUAL

IMPORTANT SAFETY INSTRUCTION

Before using the neck and shoulder massager (hereinafter - the massager), read this instruction carefully and keep it for further reference. Use the unit for intended purposes only as specified in this manual. Mishandling the unit can lead to its breakage and cause harm to the user or damage to his/her property.

- It must be provided in the transformer (DC 12V 2A) or confirm the rated voltage and then use.
- Make sure that products and transformer correctly contact.
- Do not pull the connection by force, the bad connection of transformer will cause the product to stop work.
- Do not use it in the bathroom and humid environments.
- Do not immerse the massager into water or other liquids.
- If the product falls into the water, do not touch the wet product.
- Please do not place the product on the solid surfaces to use, because it will be damaged.
- To avoid overheating, electric shock or injury, do not use this massager's close to pillows or blankets.
- Do not leave the operating unit unattended.
- Handle the massager with care, do not drop it, do not expose it to direct sunlight, high temperature and relative humidity.
- This unit is not intended for use by children.
- Do not leave children unattended to avoid them using the unit as a toy.
- The unit is not intended to be used by people with physical, sensory or mental disabilities (including children) or by persons lacking experience or knowledge if they are not under the supervision of a person who is responsible for their safety or if they are not instructed by this person on the usage of the unit.
- To avoid damages, transport the unit in the original package only.
- If a fault occurs, do not attempt to repair the massager by yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- If the power supply cable is damaged, it may only be replaced by an authorised dealer or by qualified electricians. In order to prevent hazards, always send the massager to the service centre.

SPECIFICATIONS

Rated voltage: DC 12V

Rated power: 24W

Rated frequency: 50HZ

OPERATION

- Connect the adapter jack into the adapter jack connector, ensure that the adapter jack is securely connected to the adapter jack connector.
- Connect the power adapter to the electrical outlet.
- Press the On/Off button to turn on the device.
- Apply the massager to targeted area.
- Press the Adjustable intensity control  to select the intensity.
- Press the  to select the rotation direction.
- Press the Heat button  to turn on the heat function.
- After use, turn the device off.

NAMES OF PARTS AND FUNCTIONS


— ON/OFF button

— Positive/reverse rotation

— Adjustable intensity control

— Heat function

SUGGESTED USAGE

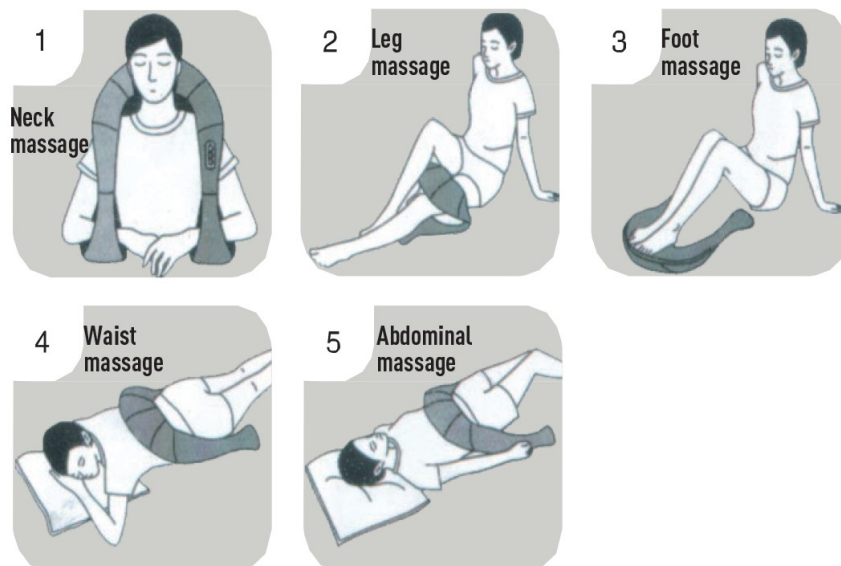CLEANING AND CARE

- Before cleaning, ensure that the main switch is turned off and the power adapter is disconnected from the electrical outlet. The main unit structure and fabric can be cleaned thoroughly with a piece of soft, dry cloth.
- Do not use water (eg. a wet cloth) to clean the massager, it may result in malfunction.
- Do not use abrasives and any solvents.

STORAGE

- After cleaning, store the massager in its original packaging and keep it in a clean, dry place.

TROUBLESHOOTING

If you encounter any problems when using this product, please refer to the following tips.

Condition	Possible reasons/Solutions
Unable to start the massager.	Make sure the power plug is fully in the socket. Confirm there is no faulty wires.
When using, it suddenly stops to work.	Check the power supply. The massager was overheating. Wait till it cools down before trying to power on again.

DISPOSAL

- This product must not be disposed of together with domestic waste.
- All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.
- Consult your local authority or your supplier for information about disposal.

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