&me.

INSTRUCTION MANUAL LEG MASSAGER

10958212

INTRODUKTION

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Vi tager forbehold for trykfejl.

INTRODUCTION

Thank you for choosing a product from &me. You have acquired a quality product which meets all applicable performance and safety standards in the EU. In order to ensure proper treatment and a long product life, we recommend that you comply with the instructions below.

Please read the manual and especially the safety instructions carefully before using the product. Keep this manual and provide it for referral along with the product when others are using it.

GENERAL SAFETY INSTRUCTIONS

- Improper use of the product can cause personal injury and damage.
- Only use the product for what it is intended. The manufacturer is not liable for damages arising from improper use or handling.
- Always place the product on a dry, flat and stable surface well away from table edge or similar.
- The product is intended for household use only. The product is not for professional use.
- The product engine/electronics, cord or plug cannot be submerged in water, and no water may get in contact with the engine/electronics.
- Never touch the product, cord or plug with wet or damp hands.
- Do not expose the product to direct sunlight, high temperatures, humidity, dust or corrosive substances.
- Never leave the product when it is turned on.
- When the product is in use, it should be kept under constant supervision. Children should always be supervised when using the product to ensure that they do not play with the product. The product is not a toy.
- The product may have sharp parts and must, therefore, be kept out of reach of children.

- Only use the accessories that come with the product or that are recommended by the manufacturer.
- The product must not be used by people with reduced sensitivity, physical or mental disability, or people who are not able to operate the product unless they are supervised or instructed by a person responsible for their safety.

CORD AND PLUG

- Do not let the cord hang over the edge.
- Keep the cord and the product away from heat sources, hot objects and open flames.
- Unplug the power cord before cleaning and if the product is not in use. Do not pull the cord when removing the plug from the socket, rather grasp the plug itself.
- Make sure that the cord is fully extended.
- Check regularly that the cord and the plug are not damaged and do not use the product if this is the case, or if it has been dropped on the floor, dropped in water or damaged in another way.
- If the product or plug is damaged, the product must be inspected and, if necessary, repaired by an authorised repairer. Never try to repair the product yourself. Please, contact the store to see if there is a repair that falls under warranty. Tampering of the product will void the warranty.

RULES OF COMPLAINTS THE WARRANTY DOES NOT APPLY:

- If the above is not observed.
- If there has been unauthorised intervention of the product.
- If the product has been mishandled, subjected to rough treatment or some other form of damage.

PRODUCT DISPOSAL

Used product should not be disposed as household waste. According to legislation, the product should be disposed of properly so that the materials can be recycled in an environmentally

friendly manner. Further information is available from the responsible local authority or local disposal company. The product is equipped with the following logo:


IMPORTANT SAFETY INSTRUCTIONS:

All safety and operating instructions must be read, understood, and followed before operating the unit.


CAUTIONS:

- 1. Seniors and/or first-time users should limit usage to 5 minutes per session. Gradually increase the session duration by 5 minutes as your body becomes accustomed to the device.
- 2. Never use the product in bathrooms or wet environments.
- 3. Store the product in a cool, dry place.
- 4. Avoid holding the product by its wire to prevent malfunction.
- 5. Do not sit on the cushion to avoid damage.
- 6. Wait at least one hour after meals before using the device.
- 7. Children or individuals with physical disabilities should only use the device under adult supervision.
- 8. Do not attempt to repair the device yourself or seek repairs from unauthorized persons or agents.

DESCRIPTION OF EACH PART

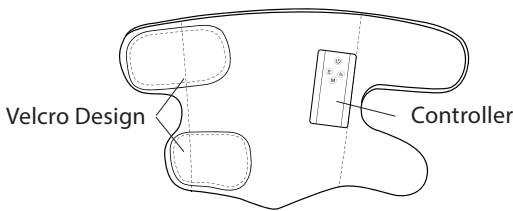CONTROLLER

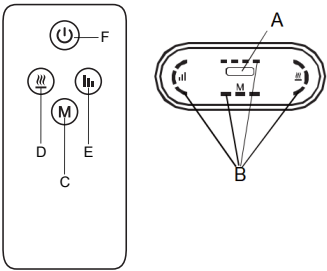- A. Type-C Charging Port
- B. Indicator Light
- C. Mode Button
- D. Heating Button
- E. Intensity Adjustment Button
- F. Power Button

USAGE

- 1. Wrap the product around your calf or arm and secure it using the Velcro straps, adjusting them to your preferred fit.
- 2. Use the controller to select your desired intensity level or mode.

CONTROLLER OPERATION

The product operates with a 3.7V/2600mAh battery.

CHARGING INSTRUCTIONS:

- During charging, the battery indicator lights will illuminate solid blue one by one.
- Once fully charged, all indicator lights will remain solid blue.

BUTTON FUNCTIONS

Power Button

- To turn on: Long press for 2 seconds. The product defaults to Mode 1 and Intensity 1, with corresponding indicator lights illuminated.
- To turn off: Long press for 2 seconds again.

MODE BUTTON

- Short press to cycle through modes:
- Mode 1 → Mode 2 → Mode 3 → Air Massage Off → back to Mode 1.

AIR INTENSITY BUTTON

- Adjusts the air intensity in three levels:
- L1 (Low): 1 indicator light on.
- L2 (Medium): 2 indicator lights on.
- L3 (High): 3 indicator lights on.
- Intensity cycles as: L1 → L2 → L3 → L1.

HEATING BUTTON

- Three heating levels to choose from:
- 50°C: 1st press.
- 55°C: 2nd press.
- 60°C: 3rd press.
- Off: 4th press.
- Heating cycles as: 50°C → 55°C → 60°C → Off.

NOTE: When the air mode is in standby, the heating function can be activated separately.

CLEANING AND MAINTENANCE

Cleaning

WARNING:

- Always switch off the device before cleaning.
- Use only a slightly damp cloth to clean the device.
- Ensure that no liquid enters the device or its accessories.
- Allow the device to dry completely before using it again.
- Do not clean the device in a dishwasher.
- Avoid using abrasive cleaning products or hard brushes.

CARE

- When not in use, store the product in a cool, dry place, preferably in its original packaging.
- Avoid placing heavy items on top of the product, as this may cause damage.

DISPOSAL

IMPORTANT:

- Do not dispose of batteries in household waste.
- Batteries may contain hazardous heavy metals and require special disposal methods.
- Please follow your local regulations for battery disposal or recycling.

TECHNICAL DATA

Article No. 10958212

Lithium Battery: 2600mAh

Input: DC5V, 2A

PRODUCER

Epiq ApS

Krøyer Kielbergs Vej 3, 3. th

8660 Skanderborg, Denmark

www.epiq.dk

We cannot be held responsible for any printing errors.